



INFRARED WAVELENGTH HELIUM-NEON LASER SYSTEM

USER'S GUIDE TO OPERATION

| UNPACKING | 3 |
|--|----|
| SHIPPING INVENTORY | 3 |
| SAFETY INFORMATION | 3 |
| Table 1: Eyewear Optical Density, Maximum Permissible Exposure, and Nominal Ocular Hazard Distance | 4 |
| COMPLIANCE | 5 |
| Figure 1: Laser Head Label Placement Diagram | 6 |
| FIGURE 2: POWER SUPPLY LABEL PLACEMENT DIAGRAM | 6 |
| OPERATING INSTRUCTIONS FOR POWER SUPPLIES | 7 |
| FIGURE 3: FRONT PANEL OF POWER SUPPLY | 7 |
| FIGURE 4: BACK PANEL OF POWER SUPPLY | |
| FIGURE 5: CONFIGURING INPUT LINE VOLTAGE SELECTOR | 8 |
| FIGURE 6: CONFIRMING SELECTED VOLTAGE | 9 |
| Figure 7: Shutter Position Matrix and Illustration | 10 |
| OPERATING INSTRUCTIONS FOR LASER HEAD WITH POWER SUPPLY | 10 |
| BASIC TROUBLESHOOTING | |
| PRODUCT WARRANTY STATEMENT | 11 |
| EU DECLARATION OF CONFORMITY | 13 |

UNPACKING

- Inspect the packing container for any damage that may have occurred during shipment.
- Contact the shipper immediately if the laser or power supply has been damaged during shipping.
- Packing materials are specifically designed to protect against shipping damage.
- Carefully remove the laser system from the packing container.
- Inspect for damage including dents, scratches or broken glass within the laser housing.
- Contact retailer immediately if the laser head or power supply appears to be damaged.

SHIPPING INVENTORY

The equipment included in this shipment will match the packing slip attached to the box. Verify that the correct product was shipped to you by matching the serial number(s) on the packing list. Notify retailer immediately if the shipment is incomplete or if an incorrect item was shipped.

Note: Please keep and store all packaging materials received with your product in case a return needs to be made.

This package contains the following:

- 1 Cylindrical Laser Head
- 1 Power supply
- 2 Keys
- 1 Remote Interlock Connector
- 1 Fuse Assembly, Fuse, and Manual Voltage Selector
- 1 115V Power Cord
- 1 230V Power Cord

SAFETY INFORMATION

The laser described in this user's guide emits infrared light. Infrared light is not visible to the human eye. This laser is safe to operate provided that the user complies with all safety warnings. It is recommended that all personnel who will operate or be in the vicinity of the laser during operation read and be familiar with this manual as well as be made aware of the following safety warnings.



- Alterations, modifications, or improper installation can result in shock, visible and/or infrared radiation exposure, and void of the warranty in part or entirety.
- Ensure that the laser head is securely connected to the power supply. To
 prevent faulty operation make sure that the male connector is fully
 seated in the back of the power supply.
- The helium neon laser can collect charge during operation acting as a capacitor. Once the power supply is turned off, best practice is to leave the high voltage cable connected to the power supply. If the high voltage

- cable is removed from the power supply the prongs should be shorted together to discharge any remaining energy. A way to accomplish this is to simultaneously touch both prongs to a conductor.
- Do not attempt to open the sealed laser housing or the power supply.
 The power supply and laser is not user accessible and service operations
 inside the enclosure must only be performed by authorized and trained
 personnel. Opening the laser or power supply will result in loss of
 warranty.
- Avoid direct exposure to the laser beam.



- This system is an infrared laser system and the emissions are not visible to the naked eye.
- Laser emissions can be hazardous to the eyes.
- Never look directly into the laser light source or at scattering laser light from reflective surfaces.
- Never sight down the beam into the source of the laser emissions.
- Install the laser so that the laser beam is not at eye level.
- Whenever the laser is operating and the beam is not in use, block the beam with the shutter on the output aperture.
- As a precaution against accidental exposure to either the laser beam or its reflection, operators need to wear laser safety glasses designed for this type of laser.

Table 1: Recommended Eyewear Optical Density, Maximum Permissible Exposure, and Nominal Ocular Hazard Distance by Part Number.

| Part Number | Description | Recommended Eyewear Optical Density (OD) | Maximum Permissible Exposure (MPE) W/cm2 | Nominal Ocular Hazard Distance (NOHD) (m) |
|---|--|--|---|--|
| 32734 | LASER SYSTEM, 633nm, STBLZD, >1.5mW | 0.28 | 0.0026 | 5.63 |
| 32172 | LASER SYSTEM, 3.39µm, 2.0mW, POL | 0 | 0.1 | 0 |
| 30990 | LASER SYSTEM, 633nm, 5.0mW, RANDOM | 0.74 | 0.0026 | 13 |
| 30025 | LASER SYSTEM, 633nm,1.5mW,POL, 9.5" | 0.3 | 0.0026 | 5.63 |
| | TUNABLE LASER, 633nm, 4.0 mW, POL | 0.74 | 0.0026 | 13 |
| | TUNABLE LASER, 612nm, 2.5 mW,POL | 0.48 | 0.0026 | 7.02 |
| 30602/30603 | TUNABLE LASER, 604nm, 0.5 mW, POL | 0 | 0.0026 | 0 |
| | TUNABLE LASER, 594nm, 0.6 mW,POL | 0.18 | 0.0026 | 3.68 |
| | TUNABLE LASER, 543nm,0.3 mW,POL | 0 | 0.0026 | 0 |
| 31007 | LASER SYSTEM, 633nm, 0.8mW, POL | 0 | 0.0026 | 0 |
| 31008 | LASER SYSTEM, 633nm, 0.5mW, RAN | 0 | 0.0026 | 0 |
| 30968 | LASER SYSTEM, 543nm, 0.5mW, POL | 0 | 0.0026 | 0 |
| 30993 | LASER SYSTEM, 633nm, 12.0mW, POL | 1.08 | 0.0026 | 19.3 |
| 33361 | LASER SYSTEM, 543nm, 1.5mW, POL | 0.3 | 0.0026 | 5.63 |
| 30967 | LASER SYSTEM, 543nm, 0.5mW, RAN | 0 | 0.0026 | 0 |
| 30988 | LASER SYSTEM, 633nm, 2.0mW, RAN | 0.48 | 0.0026 | 7.02 |
| 33141 | LASER SYSTEM, 1.52UM, 1.0mW, POL | 0 | 0.1 | 0 |
| 30992 | LASER SYSTEM, 633nm, 12.0mW, RANDOM | 1.08 | 0.0026 | 19.3 |
| 39582 | LASER SYSTEM, 594nm, 2.0mW, POL | 0.48 | 0.0026 | 7.02 |
| 39635 | LASER SYSTEM, 633nm, 17.0mW, RANDOM | 1.23 | 0.0026 | 24.8 |
| 38766 | LASER SYSTEM, 633nm, STABILIZED, AP | 0.26 | 0.0026 | 5.63 |
| 39727 | LASER SYSTEM, 633nm, STBLZD, >1.2mW | 0.15 | 0.0026 | 3.68 |
| 39581 | LASER SYSTEM, 543nm, 1.0mW, POL | 0.11 | 0.0026 | 3.68 |
| 40094 | LASER SYSTEM, 594nm, 1.0mW, POL | 0.18 | 0.0026 | 3.68 |
| 40138 | Laser System, 1.15/3.39µm, 5.0mW, POL | 0/0 | .0059/.1 | 0/0 |
| 40137 | Laser System, 1523/633nm, 1.0mW, POL | 0/0 | 0.1/.0026 | 0/3.68 |
| 40136 | Laser System, 1.15um, 1.0mW, Pol | 0 | 0.0059 | 0 |
| 40141 | Laser System, 543nm, 0.5mW, Multimode | 0 | 0.0026 | 0 |
| 14309 | LASER SYSTEM,633nm,35mW,POL | 1.58 | 0.0026 | 36.4 |
| 14354 | Laser System, 633nm, 3.0mW, ORTHOG POL, OPIS | 0.54 | 0.0026 | 9.15 |
| 39568 | Laser System, 543nm, 1.0mW, RAN | 0.11 | 0.0026 | 3.68 |
| Calculations performed by "The Evaluator" From Laser Institute of America, INC. www.lia.org | | | | |



- High voltage is present at all times when the key switch on the power supply is in the "ON" position.
- The power cord and plug are provided with a ground line. To avoid possible shock ensure that the plug is properly connected to a ground point at the electrical connection.
- Do not perform any operating or maintenance procedure that is not described in the user's manual as shock or injury may result.
- Do not operate this product if the cover has been removed.
- This product is for indoor use only. To prevent potential fire or shock hazard, do not expose the unit to any source of excessive moisture.
- Disconnect power cord before replacing fuses.
- Clean laser head and power supply with dry, soft cloth. Do not use liquids.

CAUTION: Use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure. Avoid unnecessary exposure to laser or collateral radiation that exceeds the accessible emission limits as determined by the National Center for Devices and Radiological Health (CDRH). Unauthorized modifications to the laser or power supply may cause irreversible damage to the system, as well as result in possible hazardous radiation exposure. These modifications will void all warranties.

COMPLIANCE

Each Helium-Neon laser system that is certified to be in compliance with the CDRH regulations, is equipped with a key-switch, remote interlock connector, emission indicator, time delay relay (built into the power supply), a shutter (located in the front bezel of the laser housing), and all appropriate warning labels.

Below are examples of compliance symbols and approximate placement on the laser head (Figure 1) and the power supply (Figure 2).

- A. Laser Head Manufacturing Label Includes the Following:
 - Model Number
 - Date of Manufacture
 - Minimum Specified Power Rating
 - Maximum Power Rating
 - Wavelength
 - Serial Number
 - Mandatory Compliance Symbols and Conformance Labels
- B. International Laser Beam Hazard Symbol
- C. Safety Warning and Laser Class Call-Out
- D. Laser Aperture Identification/Warning
- E. Caution Label for Voltage Input Selector
- F. Power Supply Label Includes the Following:
 - Model Number
 - Power
 - Date of Manufacture

- Serial Number
- G. Power Supply Manufacturing Label Includes the Following:
 - Model Number
 - Serial Number
 - Input Power
 - Replacement Fuse Information

Figure 1: Laser Head Label Placement Diagram


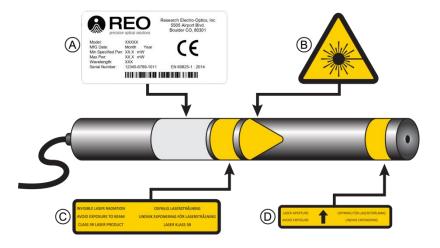


Figure 2: Power Supply Label Placement Diagram



Research Electro-Optics (REO) certifies that our laser heads and lab power supplies meet the appropriate CE requirements. For the CE regulations to be met, all lab power supplies must be used with an input line cord with a length of less than 3 meters.



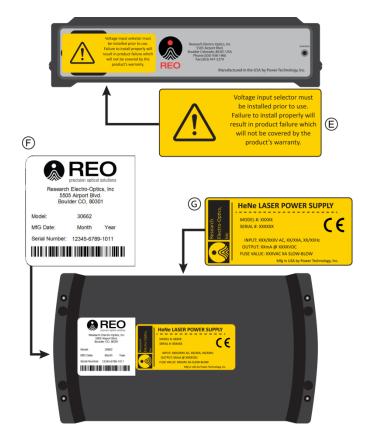

The below symbol on the product or on its packaging indicates that this product must not be disposed of with regular waste. It is the user responsibility to dispose of waste equipment according to the local laws. The separate collection and recycling of the waste equipment at the time of disposal will help to conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For further information about where the user can drop off the waste equipment for recycling, please email REO at: henelaser@reoinc.com.



OPERATING INSTRUCTIONS FOR POWER SUPPLIES

The power supply is key switch operated. It incorporates all of the appropriate safety features required by international regulations.

Figure 3: Front Panel of Power Supply

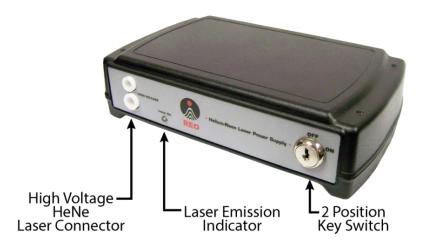


Figure 4: Back Panel of Power Supply




Fuse Assembly/Power Cord Installation

The power supply will not operate until the fuse assembly is correctly installed. The enclosed fuse assembly can be found with a 250V - 5A slow blow fuse installed. A component within this fuse assembly is the input line voltage selector. The input line voltage selector configures the correct input line voltage for the local environment. This can be accomplished by sliding the input line voltage selector into the side of the fuse assembly (Figure 5). The desired input line voltage can be seen right side up in the window on the face of the fuse assembly. Select the proper voltage for your region of the world. Once the correct line voltage is selected, insert the fuse assembly into the power cord connector located on the back of the unit. The fuse assembly will be flush with the face of the power connector when properly installed.

Note: Use only 250V-5A slow blow fuses with this product (5mm dia x 20mm length). Failing to select the correct input line voltage will result in failure of the unit and will void the warranty.

Figure 5: Configuring Input Line Voltage Selector

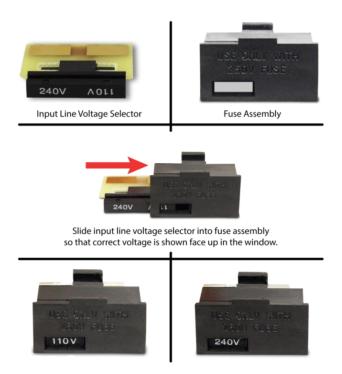

Once the fuse assembly and voltage selector are correctly installed, insert the supplied US IEC320 Female to Nema 5-15P power cord into the power supply unit. To ensure proper and optimal laser operation, only use the power supply that is appropriate for the HeNe laser head in use.

Figure 6: Confirming Selected Voltage



Fuse Replacement

To remove the fuse assembly, insert a flat blade screw driver into the groove on top of the fuse assembly and carefully remove it from the power cord connector. Remove the fuse and replace with a 250V-5A slow blow fuse (Fuse Type 2AG). While replacing the fuse assembly as described above, ensure that the input voltage selector is still properly configured.

Remote Interlock Operation

To use the remote interlock connector, remove the shorting conductor from the connector provided. The connector may then be wired into the users interlock system. The circuit is designed to activate the laser output when the two contacts of the connector are shorted together. This can be achieved via a simple switch or with relay contacts.

The voltage on this connector is nominally 5VDC and is generated internally within the power supply unit. It is completely isolated from the AC line voltage. Failure of the remote interlock feature due to improper application of external voltages is not covered by warranty.

Operating Procedures

The system is equipped with a key switch located on the front panel. This is used to activate and deactivate the laser. The emissions indicator on the front panel will come on to warn of imminent laser emission. A 10-15 second delay will begin after the key switch is turned to the "ON" position. The laser system will be "OFF" when the key is in the vertical position.


The laser emissions can be controlled from the remote interlock connector. If the factory supplied remote interlock plug is removed while the laser is in operation, laser emissions will cease immediately. If the remote interlock plug is inserted into the remote interlock connector, laser emissions will resume after a 10-15 second safety delay. If the remote interlock connector has been interfaced to external circuits, laser emissions will be controlled in a similar manner.

• The system will start within 10-15 seconds and specifications are guaranteed after thirty minutes of operation.

• Open the beam shutter located on the laser head at the output of the laser (indicated by the arrow on the aperture label).

Figure 7: Shutter Position Matrix and Illustration

| Position | Three Position Shutter | Four Position Shutter |
|----------|----------------------------|----------------------------|
| Α | Transmits 3.39µm | Transmits 6.33µm |
| В | Transmits Both Wavelengths | Transmits 1.52µm |
| С | Shut | Transmits Both Wavelengths |
| D | Not Applicable | Shut |
| | | |

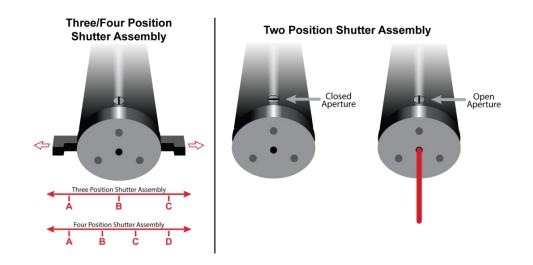

Maintenance & Service

The power supply unit contains no user serviceable parts. Return the power supply unit to retailer for any service issues that may arise.

OPERATING INSTRUCTIONS FOR LASER HEAD WITH POWER SUPPLY

The cylindrical construction of the laser head allows for easy mounting in ring clamps, V-blocks, or similar mounting hardware without affecting the laser's alignment. Use caution when clamping onto the laser as too much pressure can damage the aluminum housing or cause misalignment. The front bezel has four each 4-40 Unified National Coarse Thread (UNC) holes to secure an optional accessory mounting ring. The ring (part number 30646), which is available for purchase from Research Electro-Optics, has 1" diameter, 32 threads-per-inch (TPI) female threads.

Installation

The laser system includes a laser head and a power supply.

- Configure the AC power supply as directed in the "Operating Instructions for Power Supplies" section above.
- Insert the remote interlock connector into the back of the power supply.
- Plug the two prong high voltage cable from the laser into the front of the power supply. Confirm that the plug is well seated.
- Insert the key into the front panel of the power supply.

BASIC TROUBLESHOOTING

Once the laser and power supply have been properly connected, the laser system will operate without maintenance. Below are some guidelines if your system is not operating. When trouble-shooting take all precautions as if the laser system is fully operational.

- No emission light above the key switch after key has been turned to the on position:
 - Verify that the power cord is connected to the back of the power supply as well as a verified power source.
 - Verify that the remote interlock is fully engaged in the back panel of the power supply.
- No laser emission:
 - Visually verify that the high voltage cable is connected and well seated to the front of the power supply.
 - Verify the shutter is open.
 - o If laser does not ignite within one minute switch the control box off, wait several minutes and retry.

PRODUCT WARRANTY STATEMENT

Research Electro-Optics (REO) lasers and power supplies are warranted for twelve (12) months from the date of shipment. Items will be free from defects in material and workmanship, and will conform to the specifications provided on the REO website and published specification sheets, under normal use and service when correctly installed and maintained. REO will either repair or replace (at REO's sole option) any defective or nonconforming Product or part thereof which is returned at Buyer's expense to REO's facility, provided, that Buyer notifies REO in writing promptly after discovery of the defect or nonconformity and within the Warranty Period.

Repairs made or Products or components replaced under warranty are warranted for the remaining unexpired duration of the original Warranty Period for such Products or components. REO's standard repair charges shall be applicable forproducts returned for repair that are not covered under warranty. Out-of-warranty repairs are warranted for ninety (90) days from date of shipment of the repaired Product.

The foregoing warranty does not apply to Products or components thereof which are (a) repaired, modified or altered by any party other than REO; (b) used in conjunction with equipment not provided or authorized by REO; (c) subjected to unusual physical, thermal, electrical or optical stress, improper installation, misuse, abuse, tampering, accident, contamination, or negligence in use, storage, transportation or handling; or (d) considered a consumable item or an item requiring repair or replacement due to normal wear and tear. The foregoing warranty will also not apply if the "Warranty Void If Broken" seal located on any Product has been removed, broken or otherwise tampered with. On-site warranty repair is not covered under the foregoing warranty.

THE FOREGOING WARRANTY IS EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES. EXCEPT AS EXPRESSLY PROVIDED HEREIN, REO MAKES NO WARRANTIES, EITHER EXPRESS OR IMPLIED, EITHER IN FACT OR BY OPERATION OF LAW, STATUTORY OR OTHERWISE, REGARDING THE PRODUCTS, SOFTWARE OR SERVICES. REO EXPRESSLY DISCLAIMS ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE FOR THE PRODUCTS, OR SERVICES. THE OBLIGATIONS OF REO SET FORTH IN THIS PARAGRAPH SHALL BE REO'S SOLE LIABILITY, AND BUYER'S SOLE REMEDY, FOR BREACH OF THE FOREGOING WARRANTY. End

Warranty Procedure

Products may only be returned by Buyer when accompanied by a return material authorization number ("RMA#") issued by REO's Customer Support Center, with freight prepaid by Buyer. REO will not be responsible for any damage occurring in transit or obligated to accept Products returned for warranty repair without a RMA#. Buyer bears all risk of loss or damage to the Products until delivery at REO's designated facility. REO will pay for shipment back to Buyer for Products repaired under warranty.

Note: HeNe lasers are constructed from glass vacuum tubes and are very fragile. REO strongly recommends taking appropriate precautions in preparing your laser for shipment. Make sure the tube is well padded (equivalent to 3 inch thick high density foam on all sides) and insure that the power supply and tube are secured within the box so they do not move during shipment. REO will not assume responsibility for any shipping damages unless the laser is sent in REO-supplied packaging. REO will provide appropriate shipping materials at cost if you do not have the original shipping materials.

To request a RMA#, contact REO's Customer Support Center at:

Phone- (303) 938-1960 Fax- (303-447-3279)

Email- CustomerSupportCenter@reoinc.com

Following the receipt of your assigned RMA number return the product to Research Electro-Optics in the proper packaging with freight prepaid. Please, clearly mark the assigned RMA number on the outside of the packaging with the address as follows:

Research Electro-Optics, INCRMA# (Assigned RMA Number) 5505 Airport Blvd Boulder, CO 80301 USA



EU DECLARATION OF CONFORMITY

CE Compliance refers to following model numbers:

| Part Number | Description |
|-------------|--|
| 32734 | LASER SYSTEM, 633NM, STBLZD, >1.5MW |
| 32172 | LASER SYSTEM, 3.39UM, 2.0MW, POL |
| 30990 | LASER SYSTEM, 633NM, 5.0MW, RANDOM |
| 30025 | LASER SYSTEM, 633NM,1.5MW,POL, 9.5" |
| 31007 | LASER SYSTEM, 633NM, 0.8MW, POL |
| 31008 | LASER SYSTEM, 633NM, 0.5MW, RAN |
| 30968 | LASER SYSTEM, 543NM, 0.5MW, POL |
| 30993 | LASER SYSTEM, 633NM, 12.0MW, POL |
| 33361 | LASER SYSTEM, 543NM, 1.5MW, POL |
| 30967 | LASER SYSTEM, 543NM, 0.5MW, RAN |
| 30988 | LASER SYSTEM, 633NM, 2.0MW, RAN |
| 33141 | LASER SYSTEM, 1.52UM, 1.0MW, POL |
| 30992 | LASER SYSTEM, 633NM, 12.0MW, RANDOM |
| 39582 | LASER SYSTEM, 594NM, 2.0MW, POL |
| 39635 | LASER SYSTEM, 633NM, 17.0MW, RANDOM |
| 38766 | LASER SYSTEM, 633NM, STABILIZED, AP |
| 39727 | LASER SYSTEM, 633NM, STBLZD, >1.2MW |
| 39581 | LASER SYSTEM, 543NM, 1.0MW, POL |
| 40094 | LASER SYSTEM, 594NM, 1.0MW, POL |
| 40138 | Laser System, 1.15/3.39um, 5.0mW, POL |
| 40137 | Laser System, 1523/633nm, 1.0mW, POL |
| 40136 | Laser System, 1.15um, 1.0mW, Pol |
| 40141 | Laser System, 543nm, 0.5mW, Multimode |
| 14309 | LASER SYSTEM,633NM,35MW,POL |
| 14354 | Laser System, 633nm, 3.0mW, ORTHOG POL, OPIS |
| 39568 | Laser System, 543nm, 1.0mW, RAN |

Year C € mark verified: 2016

Type of Equipment: Electrical equipment for measurement, control and laboratory use in industrial locations.

Manufacturer: Research Electro-Optics, Inc.

5505 Airport Blvd. Boulder, CO, 80301 United States of America

Compliance was demonstrated to the following standards:

- IEC/EN 60825-1 | Laser Safety
- IEC/EN 61010-1 | Electrical Safety
- IEC/EN61326 | EMC

Note: The lasers listed in the table are designed to operate in a controlled electromagnetic environment; i.e., industrial laboratory settings. Lasers may interfere with R.F. communications such as mobile telephones and wireless networks.

Brian P. Turner Product Line Manager Research Electro-Optics, Inc. 5505 Airport Blvd Boulder, CO 80301 USA