




## **X-Cite® TURBO Field Update Instructions**

These instructions will guide you through the process of updating the settings of an X-Cite TURBO unit via a computer and USB cable.

### **Connect**

1. Connect the X-Cite unit to your PC with the supplied USB cable.
2. Turn on the X-Cite unit.
3. Let the driver install from Windows update, or download and install from:  
<https://www.excelitas.com/product/x-cite-turbo-multi-wavelength-led-illumination-system>
4. After the driver has installed, verify the connection in Device Manager.

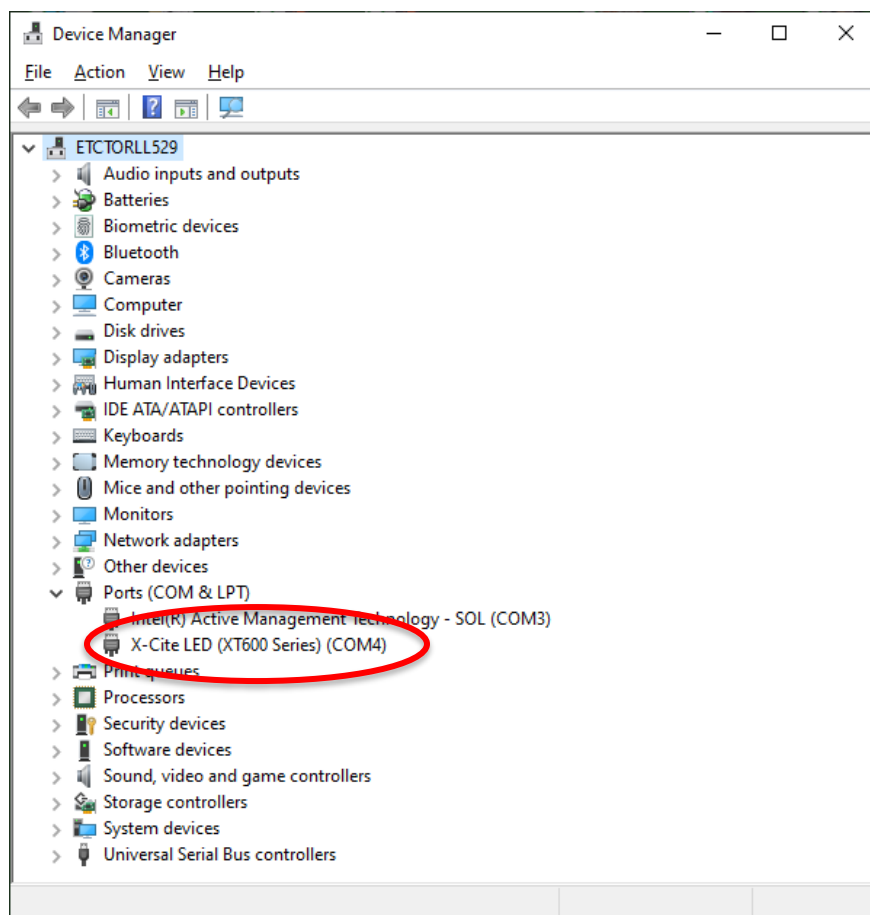


### **Obtain the Field Update Files**


5. Obtain the files from website, Tech Support or Product Management.
6. If the update zip file ends with .zi\_, save to your computer and rename .zip
7. Extract the contents of the zip file into a folder on your computer.



8. Run "Turbo-Fix.exe" from the folder created above.
9. Click "Fix" at the prompt.

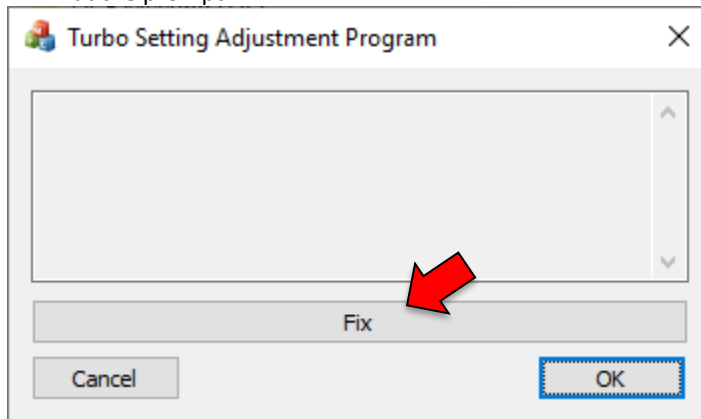

10. Click "OK" after the utility runs and dialogue box shows confirmation.

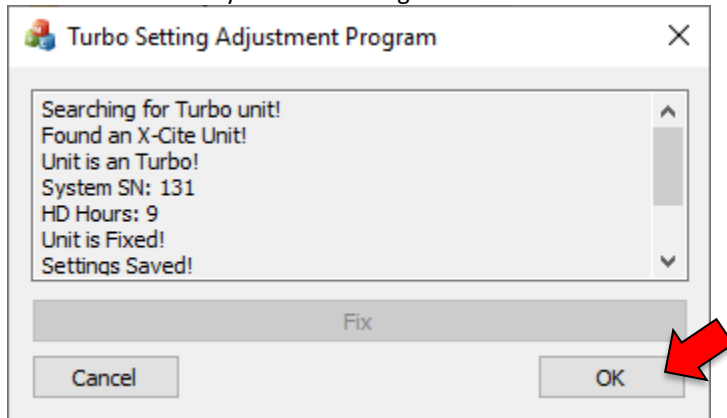

11. X-Cite TURBO is now updated with the new settings.

## Contact Information

Excelitas Technologies  
2260 Argentia Road  
Mississauga ON  
L5N6H7 CANADA

Technical Support: [techsupport@excelitas.com](mailto:techsupport@excelitas.com)  
General: [x-cite@excelitas.com](mailto:x-cite@excelitas.com)  
[www.excelitas.com](http://www.excelitas.com)